
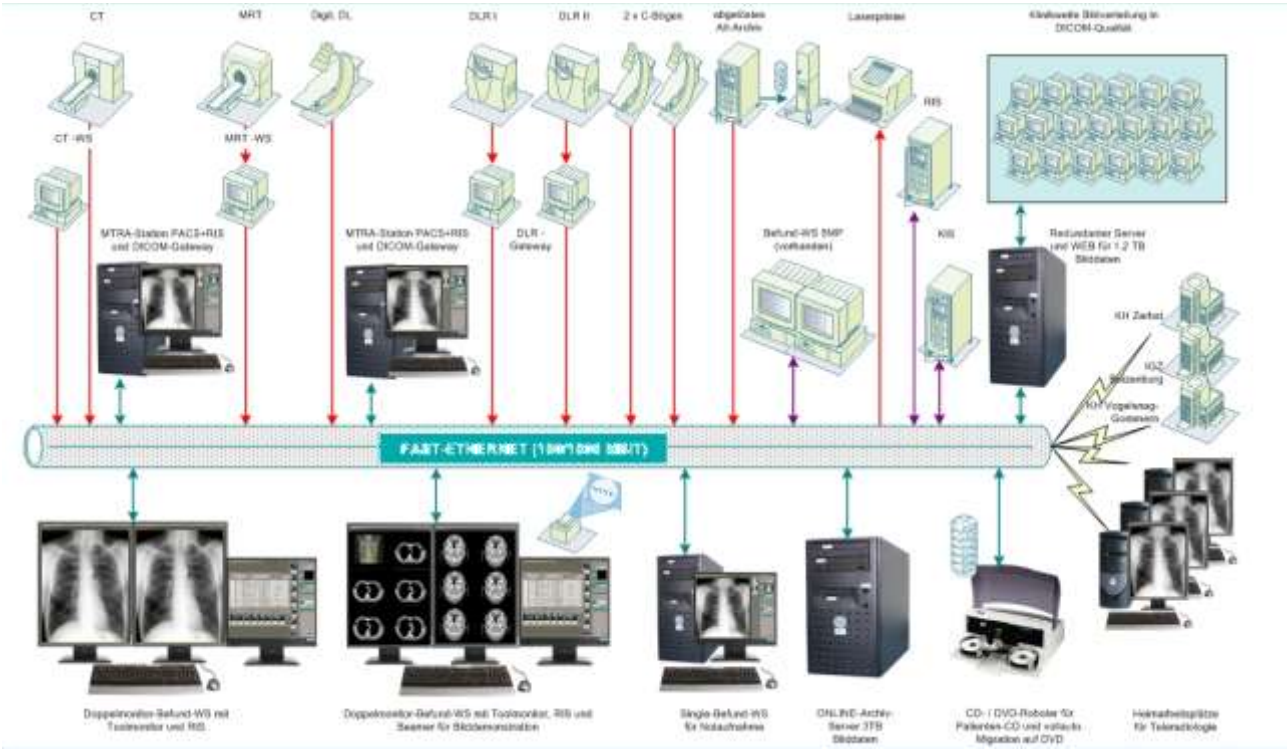


<b>Hospital/Institution Name:</b>	Kreiskrankenhaus Burg gGmbH		
<b>City/Country:</b>	Burg	Germany	
<b>Contact Persons (Phone, E-Mail):</b>	CÄ Frau Dr. Standke (Rad.)	+49 (0) 3921 96 1677	
	Herr B. Krüger (Medizintechnik)	+49 (0) 3921 96 0	
<b>Consultant (Address, Phone, E-Mail):</b>	Herr Unützer (Medbau Greifsw.)	+49 (0) 172 384 2900	
			✉ medbau@medigreif.de
<b>Connected Modalities DICOM:</b>	1 x MRT, 1 x CT, 2 x DLR, 1 x DD, 1 x WS		
<b>Connected Modalities ANALOG:</b>	2 x C-Arm		
<b>Installation date:</b>	January 2006		
<b>Investment Volume ca.:</b>	150.000 € - 180.000 €		
<b>Patients per day:</b>	<b>GB per day:</b>	120	1,1 GB
<b>Server size:</b>	<b>TB per year:</b>	1 TB	0,3 TB
<b>Number of Diagnosis Stations:</b>	8		
<b>HIS/RIS Integration with Vendor:</b>	PACS interface to SIEMENS KIS (medico//s) and RIS (Wrad)		
<b>Telecommunication with:</b>	3 x clinics and 3 x home office		
<b>Special Applications or Integrations:</b>	<p>Clinic wide image distribution via DICOM from the start</p> <p>PACS server with online RAID 1 TB – calculated for 3 years</p> <p>Redundant backup server</p> <p>All modalities with DICOM worklist – interface</p> <p>Installation time: two week during rotate clinic-work-day</p>		

## References - Kreiskrankenhaus Burg gGmbH



Krankenhaus Burg - Haus A



Krankenhaus Burg - Funktionstrakt